


How To: Use Bypass

MyChild provides a ByPass that allows you to take a patient out of a protected area without setting off alarms. You must select a departing unit and a destination unit. These can be the same. The patient status changes from “Active” to “Bypass” when a patient is in ByPass.

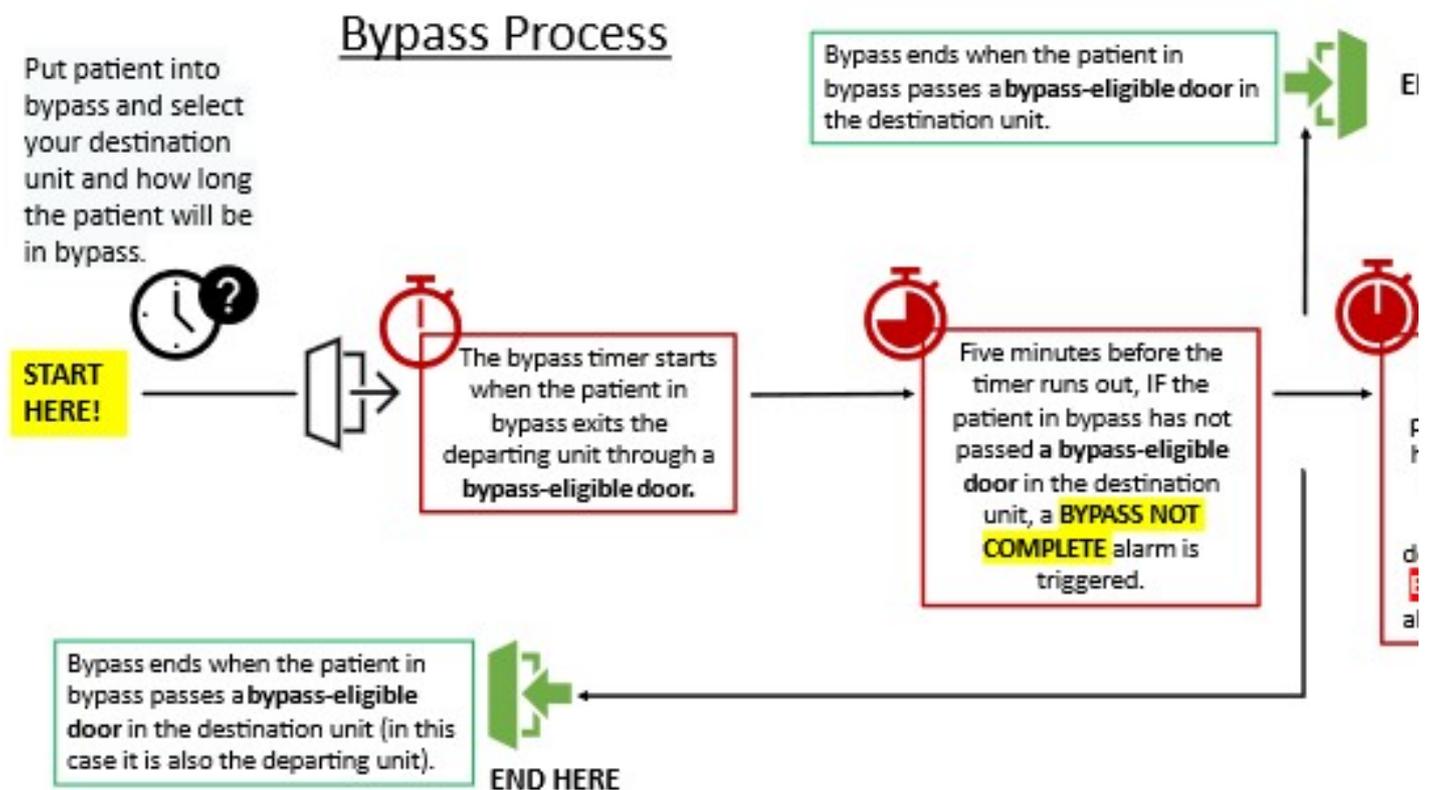

Enter a reason for the Bypass and a time limit.

The screenshot shows a 'Patient Bypass' form with the following fields:

- Departing Unit: Labor and Delivery
- Destination Unit: Mother / Baby
- Time Required (in minutes): 5
- Reason for Bypass: Moving
- Transporter Name: Jill Smith (Optional)
- Transporter Phone: (248) 555-9999 (Optional)

At the bottom of the form is a button labeled 'Put Patient In Bypass'.

Two callout boxes provide additional instructions:

- The first callout points to the 'Destination Unit' dropdown and states: "Destination unit may be the same as the departing unit OR a different protected unit."
- The second callout points to the 'Time Required' field and states: "Enter in the time you will be out of protected units." Below this text is a horizontal line, followed by: "Bypass Mode starts when the tag passes a bypass-designated door in the departing unit. Bypass Mode ends when the tag in bypass passes through a designated door in the destination unit."

Once the information is complete, click the "Put Patient In Bypass" button. This will change the patient status to Bypass. You will also see this change on the dashboard under the Patient Summary.

A level one notification will appear on the upper right corner of the screen telling you the Bypass has started once the patient leaves through the first bypass designed door.

The screenshot displays the MyChild software interface. At the top, there is a navigation bar with the following items: Home, Patients, Tags, Floorplans (with a dropdown arrow), Reports (with a dropdown arrow), Admin (with a dropdown arrow), Help (with a dropdown arrow), and Logout (with a double arrow icon). Below the navigation bar, the main heading is "LABOR AND DELIVERY".

The central part of the interface shows a floor plan of the Labor and Delivery department. The floor plan is green and white, with several blue circular markers indicating the location of patient tags. A yellow callout box is overlaid on the floor plan, containing the following information:

BYPASS STARTED


Patient Name:	Davis, Brad
Tag Number:	B2DE55
Location Start:	Labor and Delivery
Location End:	Labor and Delivery
Transporter name:	Nurse Jackie
Transporter phone #:	800-555-1234
Time of Alarm:	05/03/2022 17:43:09

At the bottom of the callout box, there are two buttons: "Accept" (green) and "Edit Bypass" (blue). Below the floor plan, there is a text box containing "Labor and Delivery".

At the bottom of the interface, there is a copyright notice: "© 2022 - McRoberts Technologies | build version: 8122".

When the time allotted for the bypass is almost expired, a "BYPASS EXPIRING" alarm is triggered.

LABOR AND DELIVERY

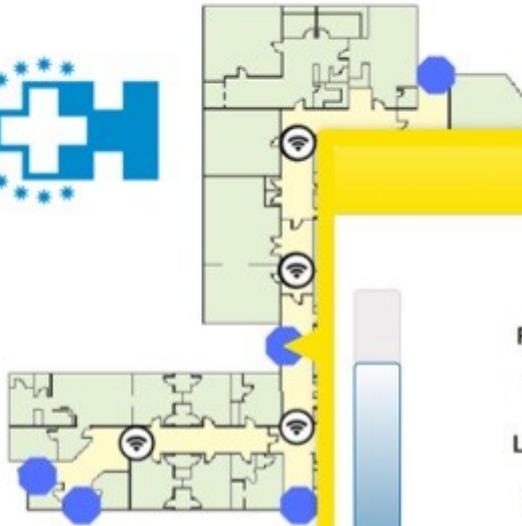

Labor and Delivery

BYPASS EXPIRING

170 secs left

Patient Name: Davis, Brad

Tag Number: B2DE55

Location Start: Labor and Delivery

Location End: Labor and Delivery

Transporter name: Nurse Jackie

Transporter phone #: 800-555-1234

Time of Alarm: 05/03/2022 17:58:06


Accept

Edit Bypass

Click "Accept" to clear the alarm from the screen. The timer will continue to run out until the bypass expires or is complete.


Click "Edit Bypass" to edit the bypass. Changes can also be made on the patient card.

When the timer reaches zero, a "BYPASS EXPIRED" alarm is triggered. All bypass alarms (except BYPASS COMPLETE) allow you to edit the bypass. Changes can also be made on the patient card.

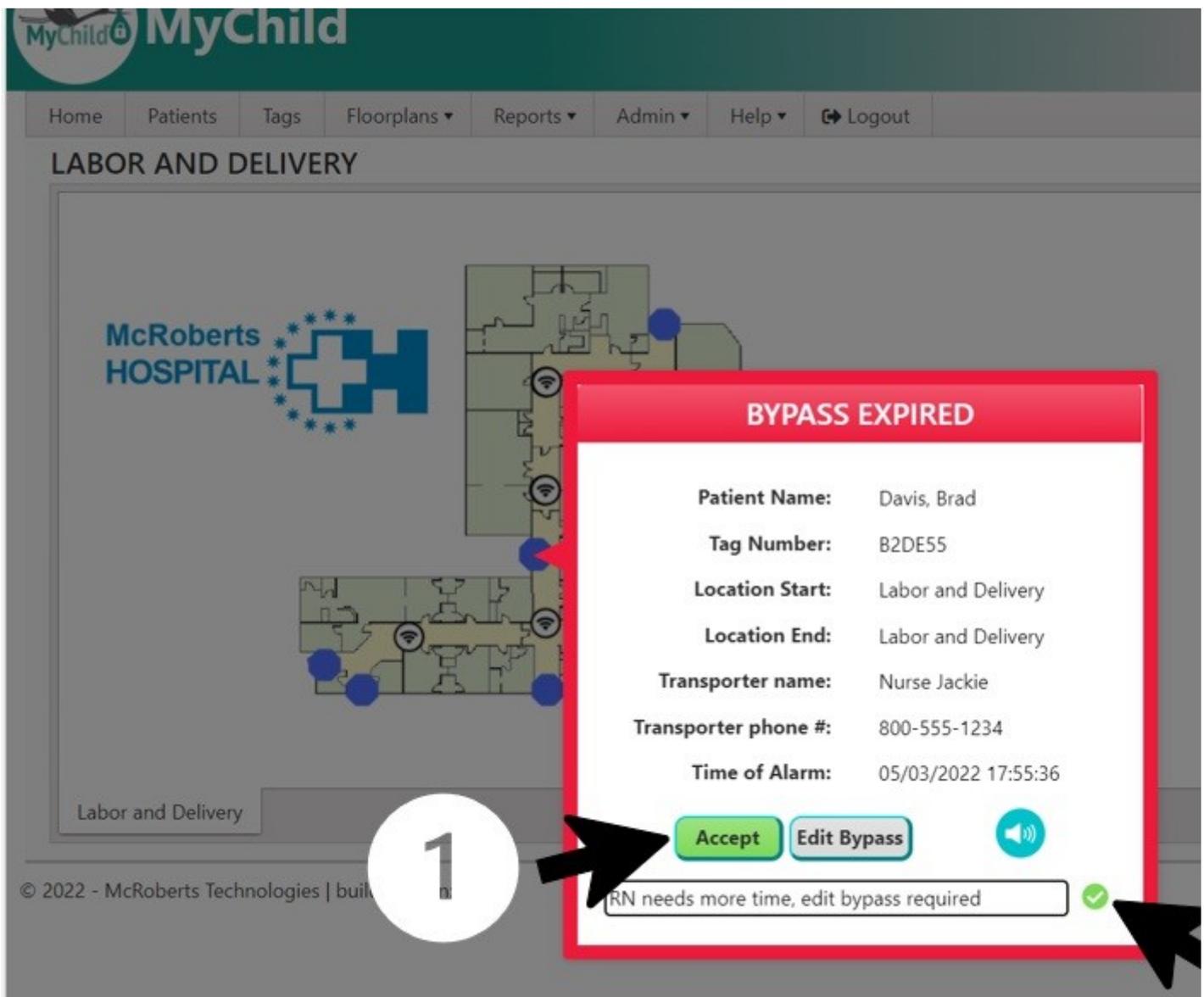

First, click the "Accept" button, then enter a message and click the green arrow to accept the alarm. The Bypass Expired alarm only happens if the time limit has run out and the patient has not passed through the second bypass designated door. If the patient arrives on time, you will receive a Bypass Complete alarm. This will show up in the upper right corner as a notification as it is a level one alarm. Click on the notification to view the alarm detail.

LABOR AND DELIVERY

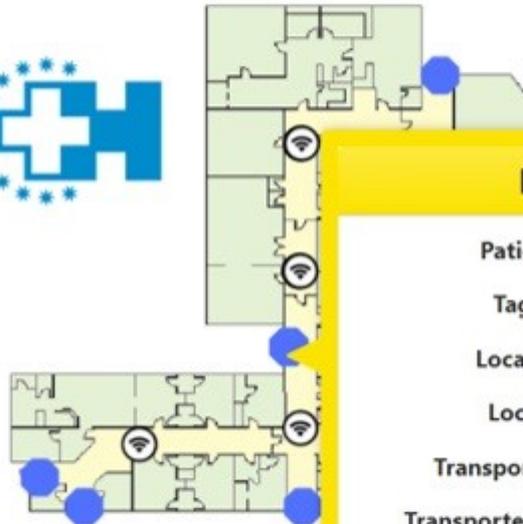

Labor and Delivery

BYPASS COMPLETE

Patient Name: Davis, Brad
Tag Number: B2DE55
Location Start: Labor and Delivery
Location End: Labor and Delivery
Transporter name: Nurse Jackie
Transporter phone #: 800-555-1234
Time of Alarm: 05/03/2022 17:44:29

When you complete a bypass by clicking “Complete,” you are acknowledging that the patient has arrived at the destination unit. The patient status will return to "Active".

For additional help, call 1-800-776-7328 and press 2 for technical support.

[Click here to download PDF.](#)

Revision #21

Created 1 February 2022 16:55:04 by Admin

Updated 2 May 2025 18:37:55 by Anna Olt