

# StaffWatch6 - Clinical

- [Staff Functions](#)
  - [How To: Enroll a Staff member](#)
  - [How To: Unenroll a Staff member](#)
  - [How To: Locate Staff](#)
- [How To Use Tags](#)
  - [How To: Check a Tag's Battery](#)
- [Alarms](#)
  - [Staff Member in Duress Alarm](#)
- [About StaffWatch](#)
  - [StaffWatch One-Sheet](#)

# Staff Functions

# How To: Enroll a Staff member

**STAFF WATCH**

Home **Staff** Help ▾ Logout

**Click Staff Tab**

**Tag Summary**

Status	Count
Assigned	9
Unassigned	15
Assigned Low-Battery	3
Unassigned Low-Battery	15
Total Low-Battery	18

To purchase tags [Click Here](#) Total Seen: 48

**Assigned Tag Count by Department**

Department	Count
Biomed	3
Emergency Department	0
Labor and Delivery	9

Total Assigned: 5

**UNACCEPTED ALARMS**

Staff Member in D

**CONTACT US**

For Tech Support, [Click Here](#)  
For Tags and Other Consumables, [Click Here](#)  
For Training, [Click Here](#)  
For Customer Service, [Click Here](#)

**Additional Offerings**

RTLS Offerings  
[Extended Warranty](#)  
[Maintenance/Prev](#)  
[Maintenance/Annu](#)  
For more informat

staging.mychildcloud.com:7270/Staff/Directory/

© 2022 - McRoberts Technologies | build version: 8144



STAFF  
WATCH

Home

Staff

Tags

Map ▾

Reports ▾

Admin ▾

Help ▾

Logout

## Staff Directory



STAFF  
WATCH

To assign a tag to  
a staff member,  
click here

Olt, Jackie

UNIT: Behavioral Unit  
TAG: A19BDC

FUNCTIONS ↻



STAFF  
WATCH

Olt, Alex

UNIT: Behavioral Unit  
TAG: A1DD03

FUNCTIONS ↻



STAFF  
WATCH

Home

Staff

Tags

Map ▾

Reports ▾

Admin ▾


Help ▾

Logout

## Staff Directory



STAFF  
WATCH

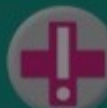

STAFF  
WATCH

Olt, Jackie

UNIT: Behavioral Unit  
TAG: A19BD

FUNCTIONS ↻

### Enroll Staff

First Name

Middle Name

Last Name

Department

Please select a unit ▾

Tag

Please select a tag ▾


Floor


Please select a floor ▾


Enroll

# How To: Unenroll a Staff member

## Staff Directory

  
NEW

  
STAFF  
WATCH

  
STAFF  
WATCH

Olt, Alex

UNIT: Behavioral Unit  
TAG: A1DD03

FUNCTIONS ↻


  
STAFF  
WATCH


Kim, Hyemi


UNIT: Behavioral Un  
TAG: A1984

Click Functions  
access all Staf  
Functions

## Staff Directory

  
NEW

  
STAFF  
WATCH

  
NO PHOTO

Olt, Alex

UNENROLL

LOCATE

  
STAFF  
WATCH

Kim, Hyemi

UNIT: Behavioral Un  
TAG: A1984



STAFF  
WATCH

Home

Staff

Tags

Map ▾

Reports ▾

Admin ▾

Help ▾

Logout

Staff Directory

+

NEW

STAFF  
WATCH

NO PHOTO

Olt , Alex

UNENROLL

LOCATE

Unenroll Staff

Are you sure you want to unenroll this staff member?

Cancel

OK

UNIT: Behavioral Unit

TAG: A1984


FUNCTIONS



# How To: Locate Staff

First, Navigate to the Staff page.

Then, click Functions to flip the card.

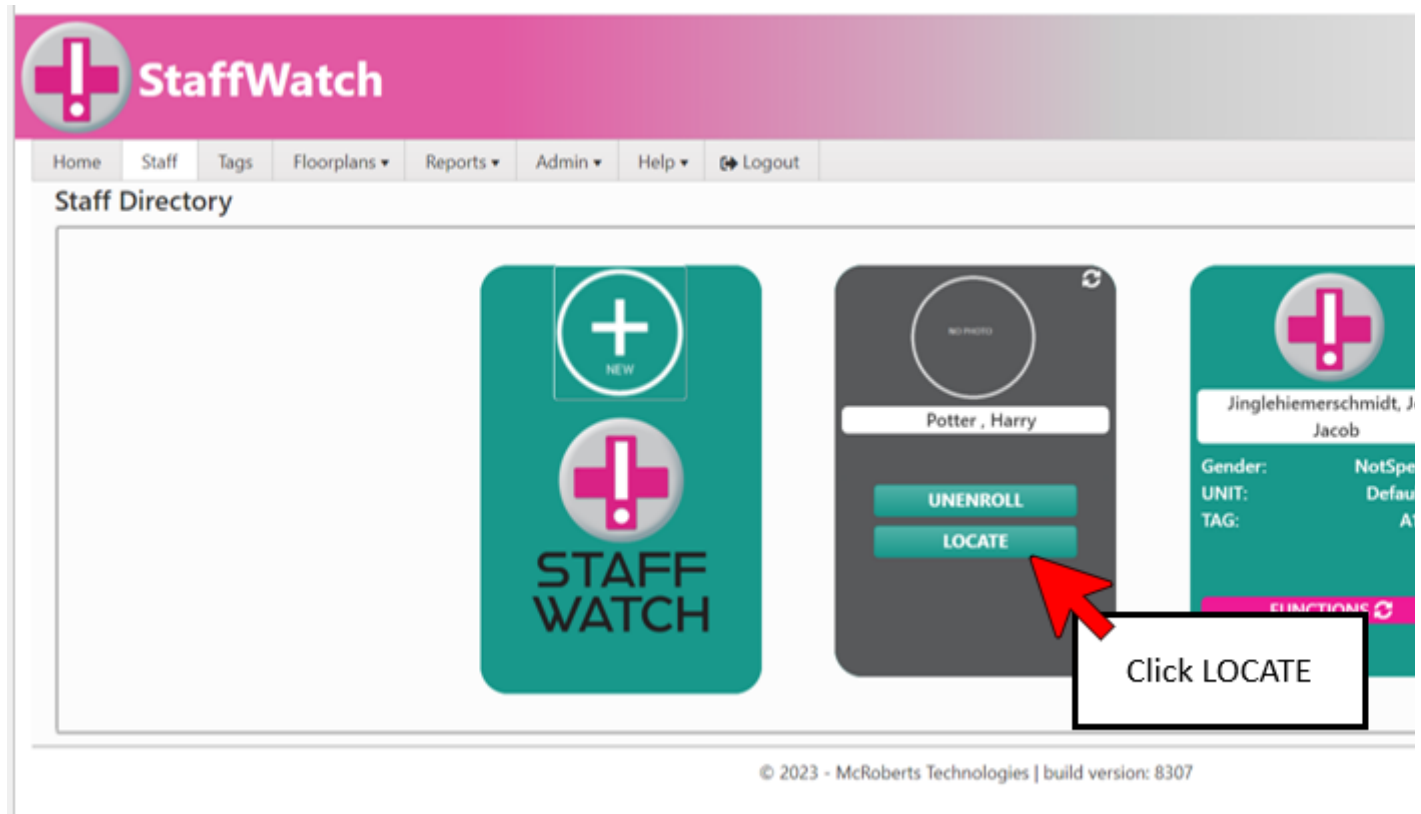

# How To Use Tags

# How To: Check a Tag's Battery

# Alarms

# Staff Member in Duress Alarm

The screenshot displays the 'STAFF WATCH' application interface. At the top, there is a navigation bar with links: Home, Staff, Tags, Map, Reports, Admin, Help, and Logout. The main content area is titled 'BEHAVIORAL UNIT' and features a map of the hospital unit. A red callout box highlights an active alarm: 'STAFF MEMBER IN DURESS'. The alarm details are as follows:

Field	Value
Device Name:	West Hall Near Rm 210
Staff Name:	Olt, Alex
Tag Number:	A1DD03
Time of Alarm:	04/20/2022 13:34:39

Below the alarm details, there is a green 'ACCEPT' button and a speaker icon. A black arrow points to the 'ACCEPT' button. A white callout box with a black border contains the text: 'Click ACCEPT to acknowledge and clear the alarm.'

At the bottom of the interface, there is a footer: '© 2022 - McRoberts Technologies | build version: 8144'.



STAFF  
WATCH

Home

Staff

Tags

Map ▾

Reports ▾


Admin ▾

Help ▾

Logout

## BEHAVIORAL UNIT

McRoberts  
HOSPITAL

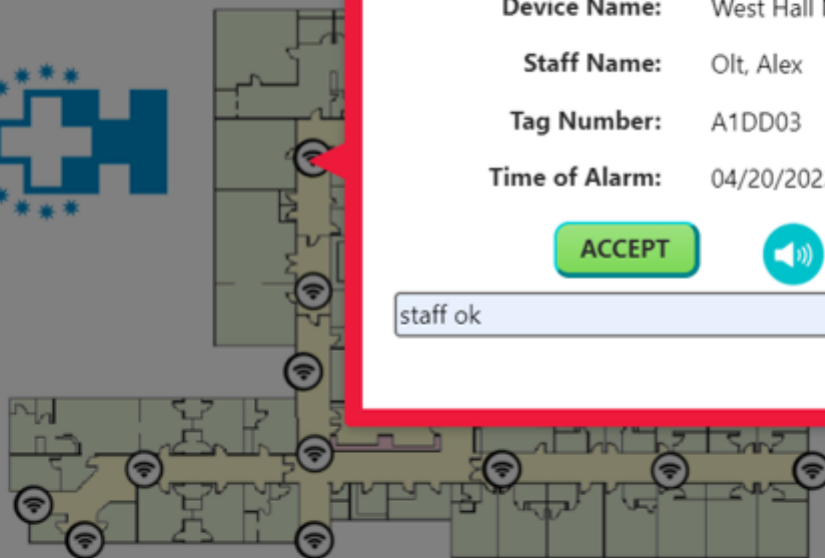

Behavioral Unit

### STAFF MEMBER IN DURESS

**Device Name:** West Hall Near Rm 210

**Staff Name:** Olt, Alex

**Tag Number:** A1DD03

**Time of Alarm:** 04/20/2022 13:34:39

ACCEPT



staff ok



# About StaffWatch

About StaffWatch

# StaffWatch One-Sheet





# STAFF WATCH


StaffWatch is an electronic panic button system for healthcare workers

FOR MORE INFORMATION CALL 1-800-776-732

OR EMAIL [SALES@MCROBERTSTECH.CO](mailto:SALES@MCROBERTSTECH.CO)

## StaffWatch® AT-A-GLANCE

- Staff members wear small, lightweight tags with panic buttons.
- Tags auto-enroll when in proximity to a receiver.
- Tag assignment to staff member is done through StaffWatch software user interface.
- When the panic button is pressed, a panic alarm is sent from the tag to the receivers via RFID and from the receivers to the server over any standard Wi-Fi network.
- Tags emit a TLM (Tag Locating Message) every 16 seconds to ensure active supervision.
- The StaffWatch server can be either cloud-based or locally-hosted.
- The panic alarm shows its location on the floor plan in the StaffWatch software application and it can be automatically sent via email or SMS to first responders.
- Locate Tag and Locate All are functions in the StaffWatch software.
- StaffWatch is cost-effective due to the low-cost, plug-in, McRoberts Wireless Receivers that create the McRoberts RTLS architecture.
- McRoberts Technologies provides installation and training and 24/7/365 technical support, including proactive remote monitoring, for life.

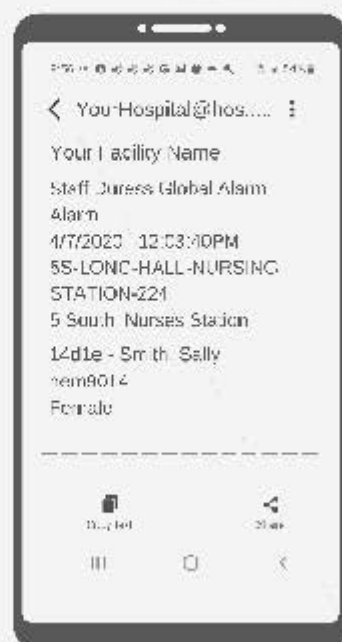

**McRoberts**  
**TECHNOLOGIES**

**WOMEN  
OWNED**

